

translumina  
therapeutics

Manufactured by :  
Translumina Therapeutics LLP  
Plot # 12 Pharma City, Selaqui, Dehradun - 248197  
t +91-135-2899795, f +91-135-2899795, www.translumina.in

Side Hole

# Trans-*Angio* TIG 2

Radial Diagnostic Catheter

**TRANSFORM** YOUR RADIAL EXPERIENCE

translumina  
therapeutics



# Trans-Angio TIG 2

## Radial Diagnostic Catheter

- Enables angiography of both RCA and LCA with one catheter<sup>1</sup> that can potentially:
  - Limit catheter exchanges
  - Shorten procedure and fluoroscopic time<sup>1</sup>
  - Lower cost per procedure
- Engages with simple clockwise and counterclockwise rotation
- Soft tip reduces potential for vessel trauma




- High Flexibility
- Superior Surface Smoothness
- Excellent Torque Control


Large Inner Lumen  
Superior Flow Rates

### Technical Data


5f	Selective
OD	0.065" / 1.65 mm
ID body	0.050" / 1.27 mm
ID tip	0.047" / 1.18 mm

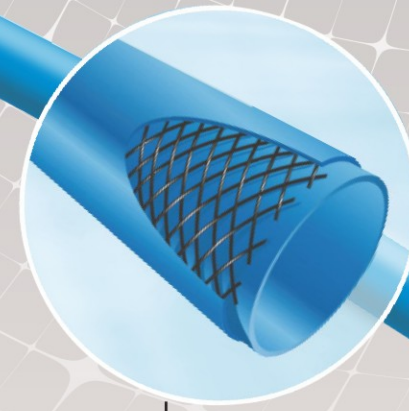

Double wire stainless steel braiding

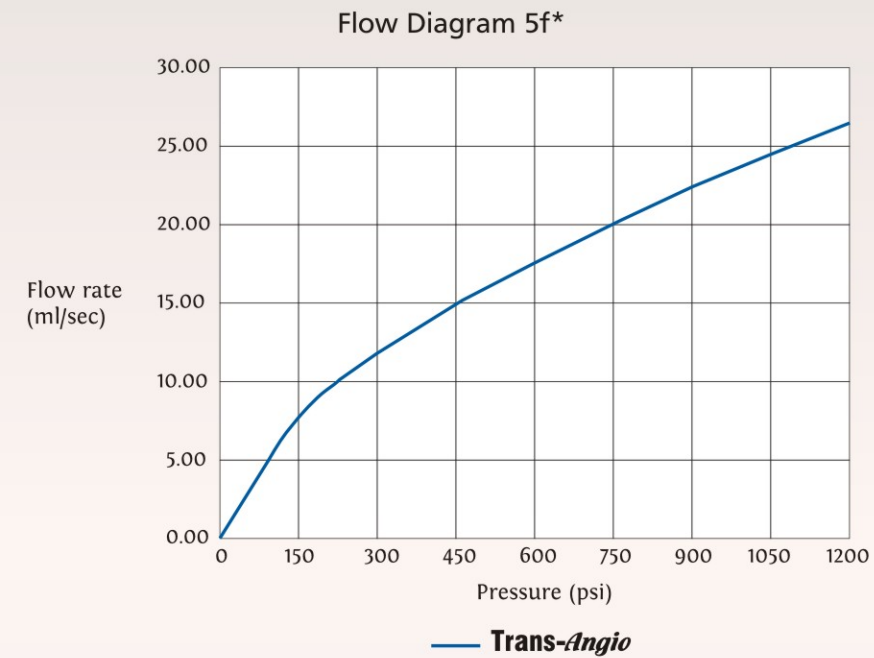

Atraumatic soft tip reduces vessel spasm

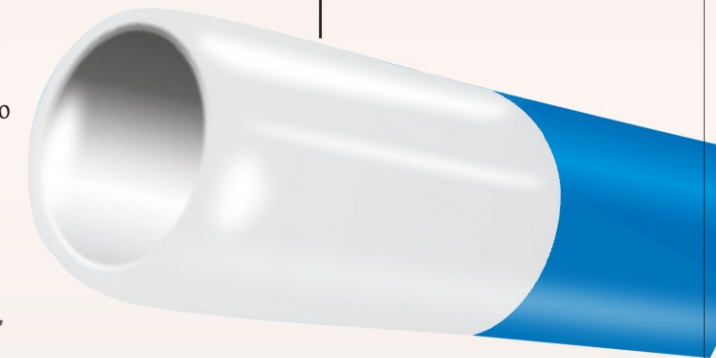

max. pressure guide wire 1200 psi / 84 bar  
0.038" / 0.97 mm

Double wire stainless steel braiding  
- Designed for 1:1 torque control